

# RC-OXYGEN NASAL CANNULA PHTHALATE FREE – INFANT



Made of soft Kraton

Package includes:

- soft Kraton nasal prongs
- tubing 2.1m




**Instructions for use –**  
Please read carefully!

**PURPOSE OF USE:** The oxygen nasal cannula is used to inhale medical oxygen through the nose. Soft Kraton nasal prongs prevent allergic reactions due to skin sensitivity or irritation at the nasal vestibule - the nose is protected. The tubing connected to the oxygen nasal cannula is kink-proof thanks to inner channels and reinforcing. Maximum flowrate 6l/min.

**USE:** Pull the oxygen nasal cannula over your head like a lariat, the nasal prongs pointing towards you. Pull the upper section forward over your ears so that the tubing is held in place by the ears like a pair of glasses. Insert the nasal prongs into each nostril. Secure the cannula in place by sliding up the fixing loop at the lower end of the oxygen nasal cannula. Connect the adaptor at the end of the tube to the oxygen source and set the prescribed oxygen flow rate.

**CLEANING AND HYGIENE:** Use a standard washing-up liquid or vinegar-and-water solution (mixing ratio 4:1) to clean the nasal prongs. Rinse with clear water. **Oxygen Nasal cannula to be used for one patient only!** Replace the oxygen nasal cannula regularly (at the latest every fortnight). There are no phthalates in the tubing and the nasal prongs.

**WARNING:** Never use nasal sprays, creams or nose drops containing oils or grease in connection with oxygen. We recommend RC-Pflege®N (Article-no. 250EN) nose and throat spray or RC-Gel NASAL (Article-no. 260). Smoking, the use of matches and open flames, such as candles or hearths, are strictly forbidden near oxygen equipment/oxygen nasal cannula.

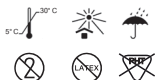

Available at pharmacies,  
health care supply stores  
or at [www.cegla-shop.com](http://www.cegla-shop.com)

**Manufacturer:**  
Respan Products Inc.  
8 Erinville Drive, Erin,  
Ontario, Canada N0B1T0

**EU representative:**  
mdi EUROPA GmbH  
Langenhagener Str. 71  
D-30855 Hannover-  
Langenhagen

**CEGLA**  
MEDICAL TECHNOLOGY

**Sale distributor:**  
R. Cegla GmbH & Co. KG  
Horresser Berg 1  
56410 Montabaur  
Germany  
Tel +49 2602 92 13-0  
[info@cegla.de](mailto:info@cegla.de)  
[www.cegla.de](http://www.cegla.de)

**REF** 712

PPN 11 03367310 92



REV.

**LOT**

**CE**  
0120